

a simple, easy, beautiful **SMILE**

the logical, innovative restorative system
designed to save time, every time



SDI


Your Smile. Our Vision.

How tooth composition determines tooth colour.

- Since the colour of the tooth is determined by the dentine, not the enamel, the dentine shade plays a major role in determining the overall tooth colour.
- The enamel is relatively translucent but requires high polishability.

Dentine is 95% of the tooth colour¹

- Single opacity
- Hue saturation
- Fluorescence

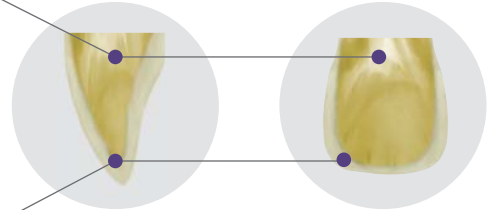

Enamel is 5% of the tooth colour¹

- Translucency
- Low hue
- Opalescence


Aura – how a linear scale provides an easy approach to tooth colour selection.

Aura is the result of five years collaboration with international materials and tooth restoration experts.


- Nearly all dentine has a similar “Hue”, falling in the red-brown-yellow range.
- Aura focuses on where tooth colours vary, the “Chroma.”
- Chroma is how strong or weak a colour (hue) is, measured in % saturation.


This linear and scientific arrangement brings predictability and simplicity to customising shades.

AURA'S EQUI-DISTANT AND BROAD CHROMA ARRANGEMENT PROVIDE PREDICTABILITY

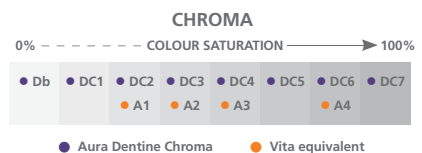

AURA IS THE ONLY SINGLE OPACITY ARRANGEMENT

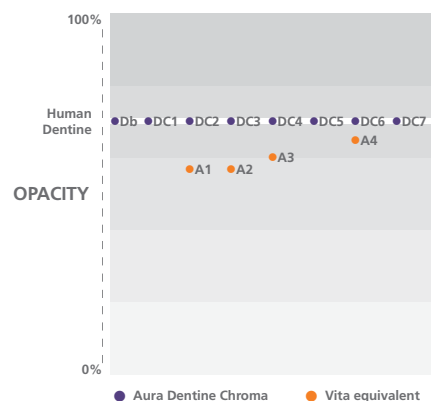

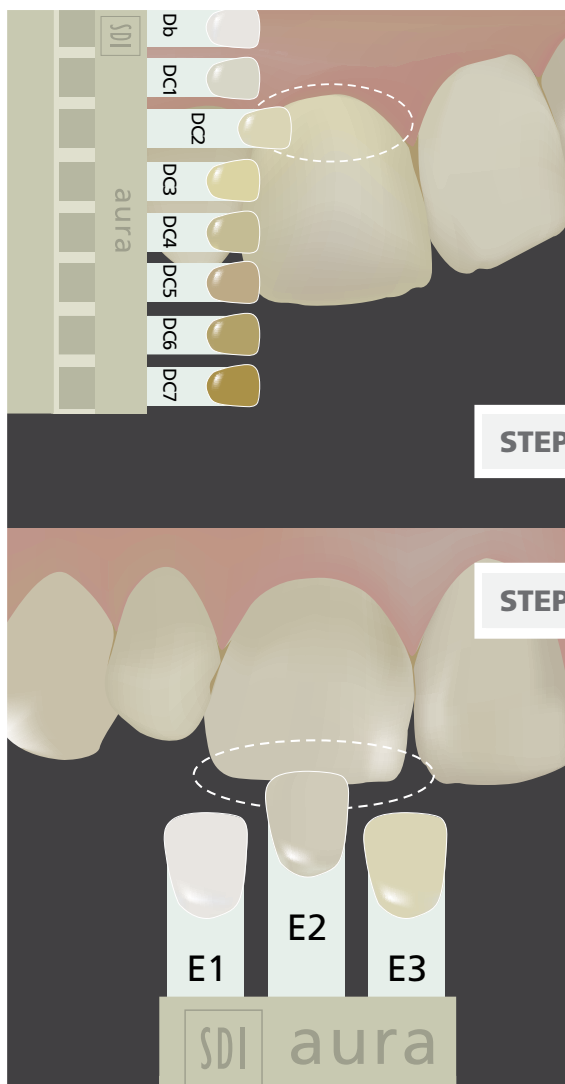

Aura – an easy approach to identifying and matching tooth colour.

Since the overall colour of the tooth is determined by the dentine, not the enamel, choosing the dentine shade should be the first priority.

STEP 1 ▶ Select the **dentine shade** where there is exposed dentine or where the cervical area of the tooth is the thinnest.

STEP 2 ▶ Select the **enamel shade** from the incisal edge of the neighbouring teeth of the tooth to be restored.

Aura – a simple approach to a mirror finish.

Aura's enamel is microfilled for easy polishability to a mirror finish. This mirror finish is retained indefinitely.

Choose the enamel shade by observing the patient:


- E1 Young or bleached teeth
- E2 Adult teeth
- E3 Senior teeth



Minimal volumetric shrinkage for maximal dimensional stability.

- Aura combines a new low-shrinkage resin composition with a new Ultra High Density (UHD) filler.


VOLUMETRIC SHRINKAGE (%)²

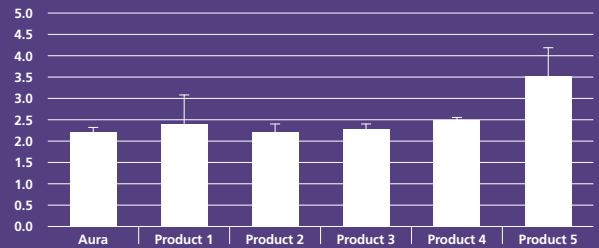

High marks for Aura's high flexural strength

- Aura's superior flexural strength means it is able to withstand forces transmitted in the mouth.


FLEXURAL STRENGTH (MPa)²

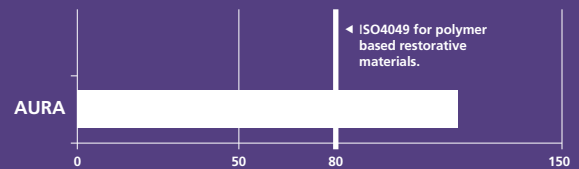

Aura's impressive compressive strength

- The unique morphology of the UHD filler provides an extremely high strength interface that can withstand high compressive forces.

COMPRESSIVE STRENGTH (MPa)²

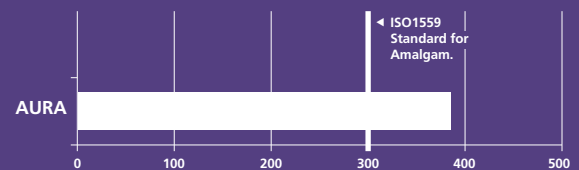

Bulk fill

- For quick posterior restorations where overlaying enamel is not necessary.
- Designed to achieve high depth of cure and strength.


Multipurpose use

- Also available in 4 multipurpose chroma shades MC2, MC3, MC4 and MC5. These are ideal for non-layering aesthetic and bulk filling requirements.

Aura, an Australian innovation.

- Based on the natural make-up and colours of enamel and dentine.
- Simple layering technique to reproduce the natural colours of teeth.
- High polishability and sculptability.
- Designed to simplify shade matching for busy clinicians.



PART NUMBER	QTY	SYRINGES	PART NUMBER	QTY	COMPLETS
 8565100		Aura Master Introductory Kit – Syringes	 8565101		Aura Master Introductory Kit – Complets
	8	4g syringe (Db, DC1, DC2, DC3, DC4, DC5, DC6, DC7)		24	0.25g complets (Db, DC1, DC2, DC3, DC4, DC5, DC6, DC7)
	3	3g syringe (E1, E2, E3)		24	0.2g complets (E1, E2, E3)
	1	4g syringe (bulk fill)		3	0.25g complets (bulk fill)
	2	AURA MASTER SHADE GUIDE		1	complet applicator
				2	AURA MASTER SHADE GUIDE
 8565102		Aura Syringe Starter Kit – Light	 8565103		Aura Complet Starter Kit – Light
	3	4g syringe (DC1, DC2, DC4)		60	0.25g complets (DC1, DC2, DC4)
	1	3g syringe (E2)		20	0.2g complets (E2)
	1	AURA LIGHT SHADE GUIDE		1	AURA LIGHT SHADE GUIDE
 8565104		Aura Syringe Starter Kit – Medium	 8565105		Aura Complet Starter Kit – Medium
	3	4g syringe (DC1, DC3, DC5)		60	0.25g complets (DC1, DC3, DC5)
	1	3g syringe (E2)		20	0.2g complets (E2)
	1	AURA MEDIUM SHADE GUIDE		1	AURA MEDIUM SHADE GUIDE
 8565106		Aura Syringe Multipurpose Kit	 8565107		Aura Complet Multipurpose Kit
	4	4g syringe (MC2, MC3, MC4, MC5)		80	0.25g complets (MC2, MC3, MC4, MC5)
	1	AURA MULTI CHROMA SHADE GUIDE		1	AURA MULTI CHROMA SHADE GUIDE
		Aura Syringe Refills			Aura Complet Refills
 8560011	1	Db – 4g	 8565011	20	Db – 0.25g
8560004	1	DC1 – 4g	 8565004	20	DC1 – 0.25g
8560005	1	DC2 – 4g	 8565005	20	DC2 – 0.25g
8560006	1	DC3 – 4g	 8565006	20	DC3 – 0.25g
8560007	1	DC4 – 4g	 8565007	20	DC4 – 0.25g
8560008	1	DC5 – 4g	 8565008	20	DC5 – 0.25g
8560009	1	DC6 – 4g	 8565009	20	DC6 – 0.25g
8560010	1	DC7 – 4g	 8565010	20	DC7 – 0.25g
8560012	1	Bulk Fill – 4g	 8565012	20	Bulk Fill – 0.25g
8560001	1	E1 – 3g	 8565001	20	E1 – 0.2g
8560002	1	E2 – 3g	 8565002	20	E2 – 0.2g
8560003	1	E3 – 3g	 8565003	20	E3 – 0.2g
8560014	1	MC2 – 4g	 8565014	20	MC2 – 0.25g
8560015	1	MC3 – 4g	 8565015	20	MC3 – 0.25g
8560016	1	MC4 – 4g	 8565016	20	MC4 – 0.25g
8560017	1	MC5 – 4g	 8565017	20	MC5 – 0.25g



Aura is proudly made in Australia by SDI Limited, Bayswater, Victoria 3153. 1 800 337 003, www.sdi.com.au
Austria 00800 022 55734, Brazil 0800 770 1735, France 00800 022 55 734, Germany 0800 100 5759
Ireland 01 886 9577, Italy 800 780625, New Zealand 0800 734 034, Spain 00800 022 55 734
United Kingdom 00800 022 55 734, USA & Canada 1 800 228 5166

References: 1. Chu S.J, Devigus A, Paravina R.D, Mieleszko A. J; Fundamentals of Color – Shade Matching and Communications in Esthetic Dentistry; p.89, 2nd Edn.
2. Data on file, SDI Limited 2014.



M100169